



LIFESENSE®

Bluetooth Thermometer

Infrared Device



FACT BOX

Product Code:
MD1877

Product Category:
Diagnostics

Contents:
1 x Non-Contact Device
2 x AAA Batteries
1 x Storage Bag

Warning:

Please keep this infrared thermometer out of children's reach. Medical assistance can't be replaced by the use of an infrared thermometer. This infrared thermometer is not waterproof. Keep it away from water.



Scan above to download our free
MEDICARE LIFESENSE APP
Measure, Track and Share your Results



The Medicare LifeSense® Infrared Bluetooth Thermometer is a cutting-edge device equipped with advanced features to streamline temperature monitoring.

With its seamless smartphone connectivity to the Medicare LifeSense® app, users can effortlessly track and manage temperature data.

Featuring a dual mode design for both forehead and ear measurements, this thermometer offers versatility and convenience.

FEATURES:

- Fast, 1 second reading
- Large easy-to-read HD LCD screen
- Child measurement mode option
- Dual mode design (forehead & ear)
- Auto shut-down after 10 seconds
- Audible and visual fever warning

INTENDED USE:

The Infrared Thermometers take human body temperature via the eardrum or forehead for both professional and home use.

OPERATING PRINCIPLE:

The infrared temperature sensor collects infrared energy emitted by the eardrum or the skin surface. After being focused by a lens, the energy is converted into a temperature reading by the thermopiles and measurement circuits.



@Medicarefirstaid



@Medicarefirstaid



@Medicare1staid



LIFESENSE®

Bluetooth Thermometer

Infrared Device




APP DETAILS

The Medicare LifeSense® App saves all your temperature measurements and creates a comprehensive chart to help build a picture of your long-term data.

Compare your previous numbers or share your information with your doctor, pharmacist or carer via email from the app.

The app also allows you to set reminders to check your temperature and you can also access all your information from the Medicare LifeSense® online portal at app.flemingmedical.ie.

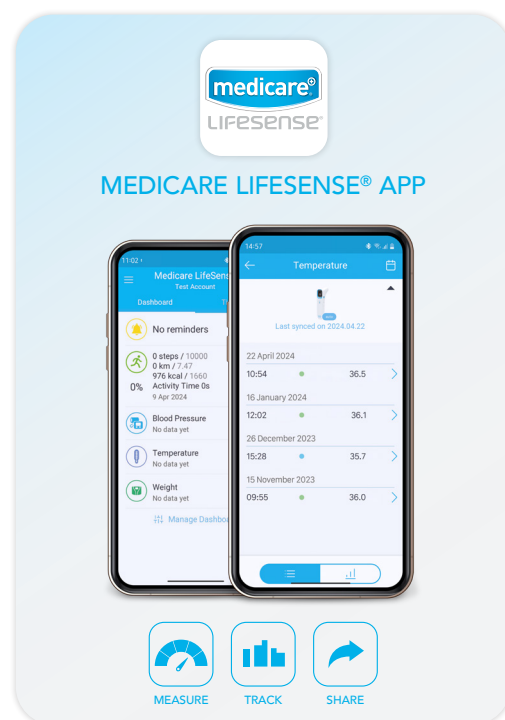

QUICK GUIDE - HEAD MODE

1



Move hair out of the way

2



Wipe away any sweat or dirt

3



Point the device at the centre of the forehead (1-3cm)



@Medicarefirstaid



@Medicarefirstaid



@Medicare1staid